

LZI Cocaine Metabolite 150 Controls

IVD For In Vitro Diagnostic Use Only



Lin-Zhi International, Inc.

| REF | | Description | Quantity |
|------|--|---|----------|
| 0347 | | Benzoylcegonine 112.5 ng/mL Level 1 Control | 1 x 5 mL |
| 0348 | | Benzoylcegonine 187.5 ng/mL Level 2 Control | 1 x 5 mL |

Intended Use

The Lin-Zhi International, Inc. (LZI) Cocaine Metabolite 150 Controls are for use as assayed quality control materials to monitor the precision of the LZI Cocaine Metabolite 150 Enzyme Immunoassay (Ref# 0340/0341) on a number of automated clinical chemistry analyzers.

Description of the Controls:

The LZI Cocaine Metabolite 150 Controls are human urine-based liquids, and ready to use. The constituent is a processed drug-free human urine matrix containing buffers, stabilizers, and less than 0.1% of sodium azide. The controls are prepared by spiking known concentrations of Benzoylcegonine into the drug-free matrix.

*Actual concentrations of these calibrators are confirmed by GC/MS or LC/MS. Values are provided only as guidelines and laboratories should determine the ranges based on their own test system and tolerance (2).

Precautions and Warning

- The LZI Cocaine Metabolite 150 Controls are for in vitro diagnostic use only. Harmful if swallowed.
- The Controls contain sodium azide, which may react with lead or copper plumbing to form potentially explosive metal azide. When disposing of such liquids always flush with a large volume of water to prevent azide build-up.
- The Controls are prepared from non-sterile human urine. They are not tested by licensed reagents for the presence of antibodies to human immunodeficiency viruses, hepatitis antigens, and/or anti-hepatitis antibodies. They should be handled as potentially infectious. Always apply good laboratory practices to avoid any skin contact or ingestion.
- Do not use the Controls beyond their expiration dates.

Preparation and Storage

The Controls are provided ready-to-use. No reconstitution is required. Label the cap before removal to identify it with the original bottle. The Controls should be stored refrigerated at 2-8°C when not in use.

Stability

When stored refrigerated at 2-8°C, the controls are stable either opened-recapped or unopened until the expiration date printed on the vial label. Store controls tightly capped when not in use. Control solution dispensed in the sample cups and left onboard the clinical analyzer should be discarded after use.

Procedure and Results

Both levels of controls (112.5 ng/mL and 187.5 ng/mL for the 150 ng/mL Cutoff) should be run daily to ensure proper assay performance. Additionally, run the controls with each new calibration and after specific maintenance or troubleshooting procedures are performed. For interpretation of results, refer to the appropriate LZI Cocaine Metabolite 150 Enzyme Immunoassay (Ref# 0340/0341) package insert (3).

Limitations

The LZI Cocaine Metabolite 150 Controls are for use with the LZI Cocaine Metabolite 150 Enzyme Immunoassay (Ref# 0340/0341) to detect benzoylcegonine in human urine only. Quality control materials should be used in accordance with local, state, and/or federal regulations or accreditation requirements.

Symbols Used

| | | | | | |
|--|------------------------------|--|------------------------------------|------|--------------------|
| | Authorized Representative | | Contents | | Reference Number |
| | Biological Risks | | Global Trade Item Number | | Safety Data Sheet |
| | Calibrator | | In Vitro Diagnostic medical device | | Temperature Limits |
| | CE Mark | | Lot Number | T.K. | Test Kit Number |
| | Consult Instructions for Use | | Manufacturer | | Use-by Date |

Bibliography

1. Urine testing for Drug of Abuse. National Institute on Drug Abuse (NIDA) Research Monograph 73, 1986.
2. Guidance for Industry and FDA Staff, Assayed and Unassayed Quality Control Material. U.S. Department of Health and Human Services. FDA, Document issued on June 7, 2007.
3. LZI Cocaine Metabolite 150 Enzyme Immunoassay (Ref# 0340/0341) package insert.

Additions, deletions, or changes are indicated by a change bar in the margin.


Notice: Adulteration of reagents, use of instruments without appropriate capabilities, or other failure to follow instructions as set forth in this labeling can affect performance characteristics, and stated or implied label claims.



Manufacturer:

Lin-Zhi International, Inc.
2945 Oakmead Village Court
Santa Clara, CA 95051
USA
Tel: (408) 970-8811
Fax: (408) 970-9030
www.lin-zhi.com

© March 2025 Rev. 3

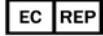

Authorized European Rep. within the EU:

CEpartner4U
Esdoornlaan 13
3951 DB Maarn
The Netherlands
www.cepartner4u.eu

