

LZI Opiate 2000 Controls



IVD For In Vitro Diagnostic Use Only



Lin-Zhi International, Inc.

REF		Description	Quantity
0337		Morphine 1500 ng/mL Level 1 Control	1 x 5 mL
0338		Morphine 2500 ng/mL Level 2 Control	1 x 5 mL

Intended Use

The Lin-Zhi International, Inc. (LZI) Opiate 2000 Controls are for use as assayed quality control materials to monitor the precision of the LZI Opiate 2000 Enzyme Immunoassay (Ref# 0330/0331) on a number of automated clinical chemistry analyzers.

Description of the Controls:

The LZI Opiate 2000 Controls are human urine-based liquids, and ready to use. The constituent is a processed drug-free human urine matrix, containing buffers, stabilizers, and less than 0.1% of sodium azide. The controls are prepared by spiking known concentrations of morphine into the drug-free matrix.

*Actual concentrations of these calibrators are confirmed by GC/MS or LC/MS. Values are provided only as guidelines and laboratories should determine the ranges based on their own test system and tolerance (2).

Precautions and Warning

- The LZI Opiate 2000 Controls are for in vitro diagnostic use only. Harmful if swallowed.
- The controls contain sodium azide, which may react with lead or copper plumbing to form potentially explosive metal azide. When disposing of such liquids always flush with a large volume of water to prevent azide build-up.
- The controls are prepared from non-sterile human urine. They are not tested by licensed reagents for the presence of antibodies to human immunodeficiency viruses, hepatitis antigens, and/or anti-hepatitis antibodies. They should be handled as potentially infectious. Always apply good laboratory practices to avoid any skin contact or ingestion.
- Do not use the controls beyond their expiration dates.

Preparation and Storage

The controls are provided ready-to-use. No reconstitution is required. Label the cap before removal to identify it with the original bottle. The controls should be stored refrigerated at 2-8°C when not in use.

Stability

When stored refrigerated at 2-8°C, the controls are stable either opened-recapped or unopened until the expiration date printed on the vial label. Store controls tightly capped when not in use. Control solution dispensed in the sample cups and left onboard the clinical analyzer should be discarded after use.

Procedure and Results

Both levels of controls (1500 ng/mL and 2500 ng/mL for the 2000 ng/mL cutoff) should be run daily to ensure proper assay performance. Additionally, run the controls with each new calibration and after specific maintenance or troubleshooting procedures are performed. For interpretation of results, refer to the appropriate LZI Opiate 2000 Enzyme Immunoassay (Ref# 0330/0331) package insert (3).

Limitations

The LZI Opiate 2000 Controls are for use with the LZI Opiate 2000 Enzyme Immunoassay (Ref# 0330/0331) to detect opiates in human urine only. Quality control materials should be used in accordance with local, state, and/or federal regulations or accreditation requirements.

Symbols Used

	Authorized Representative		Contents		Reference Number
	Biological Risks		Global Trade Item Number		Safety Data Sheet
	Calibrator		In Vitro Diagnostic medical device		Temperature Limits
	CE Mark		Lot Number	T.K.	Test Kit Number
	Consult Instructions for Use		Manufacturer		Use-by Date

Bibliography

1. Urine testing for Drug of Abuse. National Institute on Drug Abuse (NIDA) Research Monograph 73, 1986.
2. Guidance for Industry and FDA Staff, Assayed and Unassayed Quality Control Material. U.S. Department of Health and Human Services. FDA, Document issued on June 7, 2007.
3. LZI Opiate 2000 Enzyme Immunoassay (Ref# 0330/0331) package insert.

Additions, deletions, or changes are indicated by a change bar in the margin.


Notice: Adulteration of reagents, use of instruments without appropriate capabilities, or other failure to follow instructions as set forth in this labeling can affect performance characteristics, and stated or implied label claims.



Manufacturer:

Lin-Zhi International, Inc.
2945 Oakmead Village Court
Santa Clara, CA 95051
USA
Tel: (408) 970-8811
Fax: (408) 970-9030
www.lin-zhi.com

© March 2025 Rev. 3

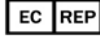

Authorized European Rep. within the EU:

CEpartner4U
Esdoornlaan 13
3951 DB Maarn
The Netherlands
www.cepartner4u.eu

